

PRODUCT DISPLAY




PRODUCT SIZE



PRODUCT DETAILS



PRODUCT PHOTOGRAPH

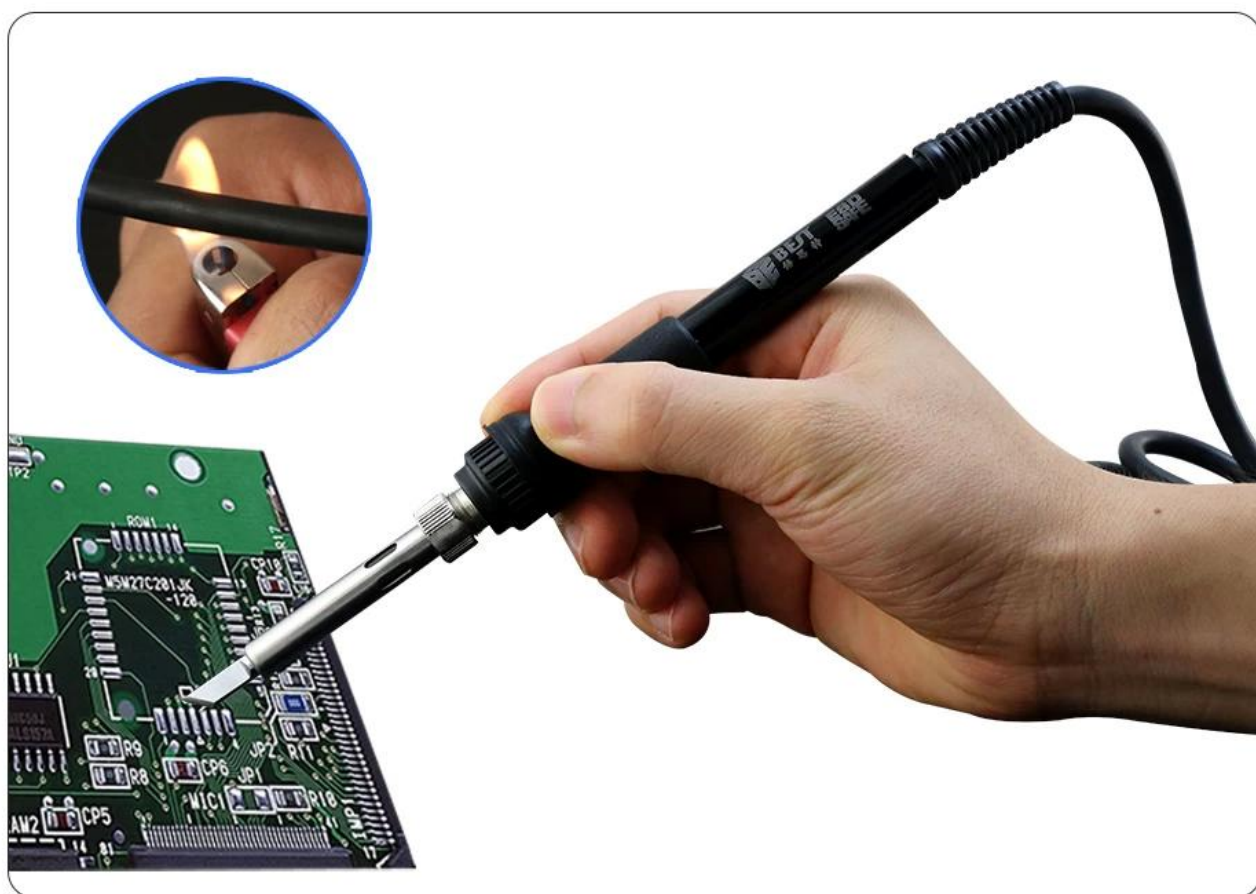

User Manual

Thank you for purchasing this precise trade from a brand of reputation. We assure you this electronic scale is a useful tool for all laboratory or business use. This electronic scale is a useful tool for all laboratory or business use. This electronic scale is a useful tool for all laboratory or business use.

I. OPERATION

- Place scale on a horizontal flat surface. Press [ON/OFF] key to turn the scale.
- Wait until "0.0" is displayed.
- Put the object on the weighing units.
- Using the [MODE] key, you can switch between the weighing units.
- Using the [COUNT] key, you can turn on or off the Backlight.

II. TARE WEIGHING

- Turn on scale as described above.
- Place the tare item on the platform.
- Press [TARE] key and "0.0" is displayed.
- Add the "net-weight" item.

III. COUNTING FUNCTION

- Turn on the scale and wait until the LCD displays "0.0".
- Place a single counting object (100g, 500g or 1000g) of the quantity you wish to weigh and count onto the platform.
- Press and HOLD [MODE] key for 3 seconds, the LCD will display a random number for example "150".
- Press [MODE] key to select the sample amount of the item that you put onto the platform.
- Press [PCS] key to confirm the above information, the LCD will display "OK" in the corner to show that the scale is in counting mode.
- Remove the sample quantity from the platform.
- You can set place any amount of your item(s) on the scale and the LCD will display a count.
- You can switch between weighing mode and counting mode by pressing [MODE] or [PCS] key repeatedly.

IV. CALIBRATION

- Insert the coin by pressing [ON/OFF] key.
- Wait until the LCD displays "0.0" and then press and HOLD [MODE] key for 3 seconds, the LCD will display "CAL" and then release the [MODE] key to repeat calibration.
- Press [MODE] key again the LCD will display "CAL" followed by the required calibration weight.
- Add the required calibration weight, the display will show "PASS" and return to normal weighing mode.
- Remove the calibration weight and press [ON/OFF] to turn off the scale and the calibration is finished.

V. SPECIFICATION:

Calibration	Auto calibration
Tare range	Tare full capacity
Auto-off	30 seconds
Operation temp.	50-85 degrees Fahrenheit
Display	7 x 3.3 digits